

Catalogue No. AB8181-200

Qty: 600 µg (3 mg/ml)

tdTomato Polyclonal Antibody

Source: Goat

General description: Goat polyclonal antibody to tdTomato (red fluorescent protein). tdTomato protein is derived from DsRed, an engineered red fluorescent protein from so-called disc corals of the genus *Discosoma*. It is a genetic fusion of two copies of the dTomato gene, which has been specifically designed for low aggregation. It's brightness and emission wavelength, makes it ideal for live animal research.

Alternative names: Cherry fluorescent protein; DsRed, mCherry, red fluorescent protein, RFP antibody.

Form: Polyclonal antibody supplied as a 200 µl (3 mg/ml) aliquot in PBS, 20% glycerol and 0.05% sodium azide. This antibody is epitope-affinity purified from goat antiserum.

Immunogen: Purified recombinant peptide produced in *E. coli*.

Specificity: In 293HEK cells transfected with cds plasmid detects a band of 55 kDa by Western blot. It also detects tdTomato in brain sections by IHC. This antibody is specific for tdTomato and mCherry proteins. It does not cross-react to GFP (green fluorescent protein).

Reactivity: Reacts against tdTomato, mCherry, red fluorescent protein (DsRed) from so-called disc corals of *Discosoma* species.

Sample	Western blot	Immuno-fluorescence	Histochemistry (paraffin)	Histochemistry (frozen)
Transfected cells	+++	+++	+++	+++

+++ excellent, ++ good, + poor, ND not determined

Usage: Western blot 1:500-1:2,000
 Immunofluorescence 1:50-1:500
 Immunohistochemistry (paraffin) 1:50-1:500
 Immunohistochemistry (frozen) 1:50-1:500

Storage: Store at -20 C for long-term storage. Store at 2-8 C for up to one month.


Special instructions: Avoid freeze/thaw cycles.

References:


1. Sada N, Fujita Y, Mizuta N, et al. *Cell Death Dis* 2020 Aug. PMID: 32811822
2. Aboualzadeh E, Phillips MJ, McGregor JE, et al. *Stem Cell Reports* 2020 Jul. PMID: 32707075
3. Laukoter S, Pauler FM, Beattie R, et al. *Neuron* 2020 Jul. PMID: 32707083
4. Muller PA, Schneeberger M, Matheis F, et al. *Nature* 2020 Jul. PMID: 32641826
5. Humbel M, Ramosaj M, Zimmer V, et al. *Gene Ther* 2020 Jul. PMID: 32632267
6. Xiong M, Tao Y, Gao Q, et al. *Cell Stem Cell* 2020 Sep. PMID: 32966778
7. Hong SJ. PhD Thesis. UC San Francisco 2020
8. Thurner M, Deutsch M, Janke K, et al. *Stem Cell Res Ther* 2020 Jun. PMID: 32532320
9. de Camargo LCB, Guaddachi F, Bergerat D, et al. *Sci Rep* 2020 Jun. PMID: 32488182
10. Morral C, Stanisavljevic J, Hernando-Momblona X, et al. *Cell Stem Cell* 2020 May. PMID:32396863
11. Edwards TJ, Fenlon LR, Dean RJ, et al. *NeuroImage* 2020 Apr. PMID:32360691
12. Nomura K, Nakanishi M, Ishidate F, et al. *Neuron* 2020 Mar. PMID:32229307
13. Zhang Y, Liu G, Guo T, et al. *Cell Rep* Mar 2020. PMID: 32234482
14. Kashiwagi M, Kanuka M, Tatsuzawa C, et al. *Curr Biol* 2020 Feb. PMID:32032507
15. Ruetten H, Wegner KA, Kennedy CL, et al. *Am J Physiol Renal Physiol* 2020 Jan. PMID:31904290
16. Best KT, Knapp E, Ketonis C, et al. *bioRxiv* 2020
17. Tran HN, Park W, Seong S, et al. *Dev Dyn* 2019 Dec. PMID:31872525
18. Stefanson O, MSc Thesis, Univ. of California - Santa Cruz, United States of America 2019
19. Pan H, Fatima M, Li A, et al. *Neuron* 2019 Sep. PMID:31324538
20. Winnubst J, Bas E, Ferreira TA, et al. *Cell* 2019 Sep. PMID:31495573
21. Park H, Choi D, Park J, et al. *Adv Sci* 2019 Sep. PMID:31763149
22. Ruud J, Alber J, Tokarska A, et al. *Cell Rep* 2019 Jun. PMID:31189104
23. Andersen RE, Hong SJ, Lim JJ, et al. *Dev Cell* 2019 May. PMID:31112699
24. Pagani M, Albisetti GW, Sivakumar N, et al. *Neuron* 2019 May. PMID:31103358
25. Ip CK, Zhang L, Farzi A, et al. *Cell Metab* 2019 Apr. PMID:31031093
26. Mathieu C, Ferren M, Jurgens E, et al. *J Virol* 2019 Apr. PMID:30728259
27. Tai Y, Gallo NB, Wang M, et al. *Neuron* 2019 Feb. PMID:30846310
28. Muller PA, Kerner Z, Pane MC, et al. *bioRxiv* 545806; Feb 2019.
29. Winnubst J, Bas E, Ferreira TA, et al. *bioRxiv* 537233; Feb 2019.
30. Albisetti GW, Pagani M, Platonova E, et al. *J Neurosci* 2019 Jan. PMID:30655357
31. Andersen RE, PhD Thesis, University of California, San Francisco, USA 2019
32. Blomfield IM, PhD Thesis, University College of London, United Kingdom 2018
33. Smidler AL, Scott SN, Mameli E, et al. *Parasit Vectors* 2018 Dec. PMID:30583744
34. Ambrozkiwicz MC, Schwark M, Kishimoto-Suga M, et al. *Neuron* 2018 Oct. PMID:30392800
35. Zalucki O, Harris L, Harvey TJ, et al. *Cereb Cortex* 2018 Oct. PMID:30272140
36. Duan L, Zhang XD, Miao WY, et al. *Neuron* 2018 Sep. PMID:30269986
37. Blomfield IM, Rocamonde B, Masdeu MDM, et al. *bioRxiv* 426015; Sept 2018
38. Niquille M, Limoni G, Markopoulos F, et al. *Elife* 2018 Mar. PMID:29557780
39. Koba S, Hanai E, Kumada N, et al. *J Physiol* 2018 Jul. PMID:30019338
40. Coux RX, Teixeira FK, Lehmann R. *Development* 2018 Apr. PMID:29511022
41. Weider M, Starost LJ, Groll K, et al. *Nat Commun* 2018 Mar. PMID:29500351
42. Unnersjö-Jess D, Scott L, Sevilla SZ, et al. *Kidney Int* 2017 Dec. PMID:29241621
43. Ishibashi F, Shimizu H, Nakata T, et al. *Stem Cell Report*. 2017 Dec. PMID:29233556
44. Lacraz GPA, Junker JP, Gladka MM, et al. *Circulation* 2017 Oct. PMID:28724751
45. Paolino A, Fenlon LR, Kozulin P, et al. *J Neurosci Methods* 2017 Sep. PMID:28917658
46. Hurni N, Kolodziejczak M, Tomasello U, et al. *Cereb Cortex* 2017 May. PMID:28334356
47. Anzo M, Sekine S, Makihara S, et al. *Genes Dev* 2017 May. PMID:28637694

References (cont.):

48. Beattie R, Postiglione MP, Burnett LE, et al. *Neuron* 2017 May. PMID:28472654

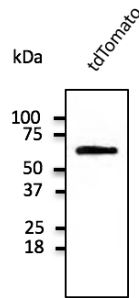

Anti-tdTomato Ab at 1/1,000 dilution; 293HEK cells transfected with cDNA; lysates at 50 µg per lane; rabbit polyclonal to goat IgG (HRP) at 1/10,000 dilution;

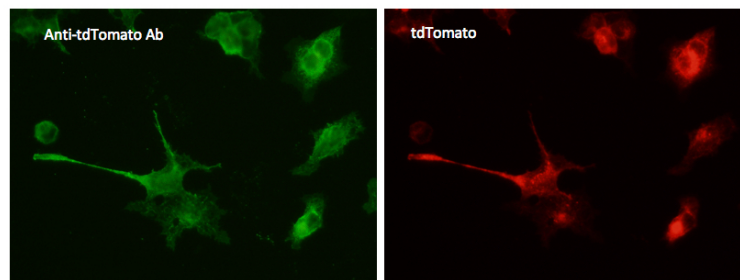

Immunofluorescence – anti-tdTomato Ab in 293HEK cells transfected with cDNA at 1/50 dilution; cells were fixed with 4% of PFA;

For research use only, not for diagnostic use

SICGEN's Proprietary Immunogen Policy

In order to produce high specific antibodies SICGEN has invested a lot of time and effort into selecting immunogen sequences. SICGEN has decided to protect this information by not publishing it on the website. However, these sequences are available on request